


# 2014 Version: Jefferson Lab Hall B Standard Equipment Manual

*Info Level 4*

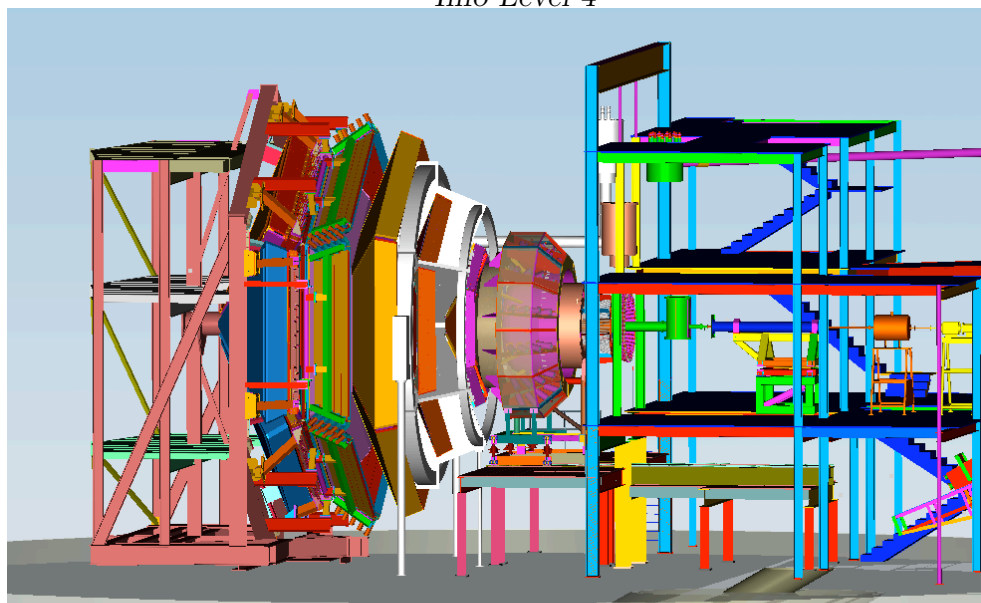

*The Hall B Collaboration*

*Editor: H. Avakian and V. Burkert*<sup>1</sup>

July 18, 2014

<sup>1</sup>Thomas Jefferson National Accelerator Facility

# Contents

<b>1</b>	<b>Hall Specific Equipment</b>	<b>4</b>
1.1	Overview . . . . .	4
1.2	Checking Tie-in To Machine Fast Shutdown System . . . . .	4
<b>2</b>	<b>General Description</b>	<b>5</b>
2.1	Introduction . . . . .	5
<b>I</b>	<b>Beamline</b>	<b>6</b>
2.1.1	Hazards . . . . .	7
2.1.2	Mitigations . . . . .	7
2.1.3	Responsible Personnel . . . . .	8
<b>II</b>	<b>HPS Spectrometer</b>	<b>9</b>
<b>3</b>	<b>HPS Chicane Magnets</b>	<b>10</b>
3.0.4	Hazards . . . . .	10
3.0.5	Mitigations . . . . .	10
3.0.6	Responsible Personnel . . . . .	10
<b>4</b>	<b>Target System</b>	<b>12</b>
4.0.7	Hazards . . . . .	12
4.0.8	Mitigations . . . . .	12
4.0.9	Responsible Personnel . . . . .	12
<b>5</b>	<b>Hall B Vacuum System</b>	<b>14</b>
5.0.10	Hazards . . . . .	14
5.0.11	Mitigations . . . . .	14
5.0.12	Responsible Personnel . . . . .	14
5.1	Silicon Tracker . . . . .	15
5.1.1	Hazards . . . . .	15
5.1.2	Mitigations . . . . .	15
5.1.3	Responsible Personnel . . . . .	16

5.2	Electromagnetic Calorimeter . . . . .	16
5.2.1	Hazards . . . . .	17
5.2.2	Mitigations . . . . .	17
5.2.3	Responsible Personnel . . . . .	17
5.3	CLAS12 Construction work . . . . .	17
5.3.1	Hazards . . . . .	17
5.3.2	Mitigations . . . . .	17
5.3.3	Responsible Personnel . . . . .	18
<b>III</b>	<b>Slow Controls</b>	<b>19</b>
<b>6</b>	<b>Overview</b>	<b>20</b>
6.1	System's Components . . . . .	20
6.2	Operating Procedures . . . . .	20
6.3	AlarmHandler . . . . .	20
6.4	StripTool . . . . .	21
6.5	Snapshot . . . . .	21
6.6	Troubleshooting Procedures . . . . .	21
6.6.1	Authorized Personnel . . . . .	21
<b>IV</b>	<b>Data Acquisition and Trigger</b>	<b>22</b>
<b>7</b>	<b>Spectrometer Data Acquisition</b>	<b>23</b>
7.1	General Computer Information . . . . .	23
7.2	DAQ checklist . . . . .	23
7.3	Beginning of Experiment Checkout . . . . .	24
7.4	Running CODA . . . . .	24
7.4.1	Some Frequently Asked Questions about DAQ . . . . .	24
7.4.2	Quick Resets . . . . .	25
7.4.3	Cold Start of CODA . . . . .	25
7.4.4	Recovering from a Reboot of Workstation . . . . .	25
7.5	Electronic Logbook and Beam Accounting . . . . .	26
7.6	Authorized Personnel . . . . .	26
<b>8</b>	<b>Trigger Hardware and Software</b>	<b>27</b>
8.1	Overview . . . . .	27
8.2	Components . . . . .	27
8.2.1	Authorized Personnel . . . . .	27
<b>9</b>	<b>Online Analysis, Data Checks</b>	<b>29</b>
9.0.2	Scaler Display and Scaler Events . . . . .	29
9.1	Analysis using CLAS12 Monitoring . . . . .	29

<i>CONTENTS</i>	3
9.1.1 Responsible Personnel . . . . .	29
<b>V Offline Analysis Software</b>	<b>30</b>
<b>10 CLAS12 Offline Analysis</b>	<b>31</b>

# Chapter 1

## Hall Specific Equipment

### 1.1 Overview

The following Hall B subsystems are considered part of the experimental endstation equipment for running the Heavy Photon Search (HPS) experiment engineering run. Many of these subsystems impose similar hazards, such as those induced by magnets and magnet power supplies, high voltage systems and vacuum systems. Note that a specific sub-system may have many unique hazards associated with it. For each major system, the hazards, mitigations, and responsible personnel are noted.

The material in this chapter is a subset of the material in the full Hall B operations manual and is only intended to familiarize people with the hazards and responsible personnel for these systems. It in no way should be taken as sufficient information to use or operate this equipment.

### 1.2 Checking Tie-in To Machine Fast Shutdown System

In order to make sure that hall equipment that should be tied into the machine fast shutdown (FSD) system has been properly checked, the hall work coordinator must be notified by e-mail prior to the end of each installation period by the system owner that the checks have been performed in conjunction with accelerator operations (i.e. checking that equipment's signals will in fact cause an FSD). These notifications will be noted in the work coordinator's final check-list as having been done. System owners are responsible for notifying the work coordinator that their system has an FSD tie-in so it can be added to the check-list.

# Chapter 2

## General Description <sup>1</sup>

### 2.1 Introduction

---

<sup>1</sup>Authors: S.Stepanyan [stepanya@jlab.org](mailto:stepanya@jlab.org)

# Part I

## Beamline

The control and measurement equipment along the Hall B beamline consists of various elements necessary to transport beam with required specifications onto the production target and the beam dump, and simultaneously measure the properties of the beam relevant to the successful implementation of the physics program in Hall B.

The beamline in the Hall provides the interface between the CEBAF accelerator and the experimental hall. All work on the beamline must be coordinated with both physics division and accelerator division in order to ensure safe and reliable transport of the electron beam to the dump.

### 2.1.1 Hazards

Along the beamline various hazards can be found. These include radiation areas, vacuum windows, high voltage, and magnetic fields.

### 2.1.2 Mitigations

All magnets (dipoles, quadrupoles, sextuples, beam correctors) and beam diagnostic devices (BPMs, scanners, beam loss monitors, viewers) necessary to transport and monitor the beam are controlled by Machine Control Center (MCC) through EPICS [1], except for specific elements which are addressed in the subsequent sections. The detailed safety operational procedures for the Hall B beamline should be essentially the same as the one for the CEBAF machine and beamline.

Personnel who need to work near or around the beamline should keep in mind the potential hazards:

- Radiation "Hot Spots" - marked by ARM of RadCon personnel,
- Vacuum in beamline tubes and other vessels,
- Thin windowed vacuum enclosures (e.g. the scattering chamber),
- Electric power hazards in the vicinity of magnets, and
- Conventional hazards (fall hazard, crane hazard, etc.).

These hazards are noted by signs and the most hazardous areas along the beamline are roped off to restrict access when operational (e.g. around the HPS chicane magnets). Signs are posted by RadCon for any hot spots. Survey of the beamline and around it will be performed before work is done on the beam line or around. The connection of leads to magnets have plastic covers for electrical safety. Any work around the magnets will require de-energizing the magnets. Energized magnets are noted by read flashing beacons. Any work on the magnets requires the "Lock and Tag" procedures [2].

Additional safety information is available in the following documents:

- EH&S Manual [2];
- PSS Description Document [3]
- Accelerator Operations Directive [4];



### 2.1.3 Responsible Personnel

The beamline requires both accelerator and physics personnel to maintain and operate. It is very important that both groups stay in contact with each other to coordinate any work on the Hall B beamline. The authorized personnel is shown in table 2.1.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
F.-X. Girod	Hall-B	6002	757-438-4523+	<a href="mailto:fxgirod@jlab.org">fxgirod@jlab.org</a>	<i>1st Contact</i>
Stepan Stepanyan	Hall-B	7196		<a href="mailto:stepanya@jlab.org">stepanya@jlab.org</a>	<i>2nd Contact</i>
Michael Tiefenback	Accel.	7430		<a href="mailto:tiefen@jlab.org">tiefen@jlab.org</a>	<i>Contact to Hall-B</i>
Hari Areti	Accel.	7187	584-	<a href="mailto:areti@jlab.org">areti@jlab.org</a>	<i>Contact to Physics</i>

Table 2.1: Hall B beamline: authorized personnel.

# Part II

## HPS Spectrometer

# Chapter 3

## HPS Chicane Magnets <sup>1</sup>

The HPS experiment will use a three magnet chicane installed on the beamline in the Hall B downstream alcove. The chicane includes the Hall B pair spectrometer magnet, an 18D36 dipole, with its vacuum chamber, and two identical H-dipoles known as "Frascati" magnets. Electron beam will pass through the magnets in the vacuum. The pair spectrometer magnet will serve as a spectrometer magnet for HPS. The two "Frascati" dipoles are to keep the deflected electron beam on the beamline to the dump. The spectrometer magnet will be powered from the Hall B pair spectrometer magnet power supply, DANFYSIK-8000. The "Frascati" magnets will be powered from so called "mini-torus" power supply, Dynapower PD42-04000103-GKLX-PY57.

### 3.0.4 Hazards

Dipole magnets can present magnetic and electrical hazards when they are energized. There is also a possible hazard of unwanted beam motion if there is a magnet power trip.

### 3.0.5 Mitigations

There are plastic covers on the connection panels for power leads on the magnets for electrical safety. Any work around the magnets will require de-energizing the magnets. Energized magnets are noted by red flashing beacons. Any work on the magnets requires a "Lock and Tag" procedure [?]. There will be beacons installed to notify when magnetic field is present. The magnet power supplies will be interlocked to the beam shutdown system (FSD). If any of power supplies will trip, beam delivery to Hall B will be interrupted.

### 3.0.6 Responsible Personnel

The chicane magnets will be maintained by the Hall B engineering group. The authorized personnel is shown in table 3.1.

---

<sup>1</sup>Authors: S. Stepanyan [stepanya@jlab.org](mailto:stepanya@jlab.org)

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Tech-on-Call	Hall-B				<i>1st Contact</i>
Doug Tilles	Hall-B	7566		<a href="mailto:tilles@jlab.org">tilles@jlab.org</a>	<i>2nd Contact</i>

Table 3.1: Hall B beamline: authorized personnel.

# Chapter 4

## Target System

The HPS target system consists of several solid material foils mounted on a target ladder. The bottom edge of the foils mounted on the ladder is free-standing so there is no thick support frame to trip the beam when the target is inserted. Target position is remotely adjustable vertically allowing different targets to be inserted. The support frame on the beam-right side of the target is made thin enough to prevent excessive flux of secondaries and radiation damage to the silicon detector in the event of an errant beam, caused, for example, by an chicane magnet trip which will move beam to the right.

### 4.0.7 Hazards

There are hazards related to moving the target frame into the beam or overheating the target foils. The stepping motor linear actuator will be operated using EPICS controls. The GUI for operation of the target will have preset coordinates for each target foil. The tungsten targets are intended to operate with beam currents up to 500 nA, which produce strong local heating. The strength of tungsten drops by an order of magnitude with temperature increases in the range of 1000 C. In addition, the material re-crystallizes above this range, which increases the tendency for cracking where thermal expansion has caused temporary dimpling.

### 4.0.8 Mitigations

There will be limit switches (hard stops) that will prevent the motion of the target ladder outside of allowed range if EPICS set values are wrong. To keep the temperature rise less than about 1000 degrees, we adjust the optics to produce an adequately large beam spot and limit the maximum current incident on the target. There will be overall beam current limit of 500 nA for the experiment.

### 4.0.9 Responsible Personnel

The target system will be maintained by the Hall B engineering group.

Name	Dept.	Phone	email	Comments
Tech-on-call	Hall-B			1st contact
D. Tilles	Hall-B		tilles@jlab.org	2nd contact
C. Field	SLAC		sargon@slac.stanford.edu	contact

Table 4.1: Personnel responsible for the target.

# Chapter 5

## Hall B Vacuum System

<sup>1</sup>

The Hall B vacuum system consists of three segments, all interconnected. The beam transport line consisting of 1.5 to 2.5 inch beam pipes, the Hall B tagger magnet vacuum chamber, and the set of vacuum chambers through the HPS detector system. Only the tagger vacuum chamber has a large window, 8 inches over 30 ft Kevlar-Mylar composite window. There is a 2.5 in diameter 7 mil Kapton window at the end of the HPS detector vacuum system, before the shielding wall, that is normally inaccessible. The vacuum in the system is provided by a set of rough, turbo, and ion pumps and it is maintained at the level of better than  $10^{-5}$  Torr.

### 5.0.10 Hazards

Hazards associated with the vacuum system are due to rapid decompression in case of a window failure. Loud noise can cause hearing loss. Also, there is a hazard related to SVT coolant leaking into analyzing magnet vacuum chamber that will degrade the vacuum and may damage readout electronics if the leak is extensive.

### 5.0.11 Mitigations

All personal working in the vicinity of the tagger vacuum chamber window are required to wear ear protection. Warning signs must be posted in that area. To mitigate a possible coolant leakage, the cooling system will be interlocked with the beamline vacuum system and the cooling system pressure gage. In an event of a leak, evident by increased vacuum system pressure or decreased cooling system pressure, the cooling system will be shutdown and valved off.

### 5.0.12 Responsible Personnel

---

<sup>1</sup>Authors: S. Stepanyan [stepanya@jlab.org](mailto:stepanya@jlab.org)

The vacuum system will be maintained by the Hall B engineering group. The authorized personnel is shown in table 5.1.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Tech-on-Call	Hall-B				<i>1st Contact</i>
Doug Tilles	Hall-B	7566		<a href="mailto:tilles@jlab.org">tilles@jlab.org</a>	<i>2nd Contact</i>

Table 5.1: Hall B beamline: authorized personnel.

## 5.1 Silicon Tracker

The silicon vertex tracker (SVT) is a compact six layer tracking system, less than a meter long, which uses silicon microstrip detectors to measure charged particle momentum and decay vertex positions. Each layer, top or bottom, consists of axial and stereo silicon sensors. The SVT is divided into top and bottom sections to avoid direct interactions with the beam and degraded electrons, and it resides in vacuum, to eliminate beam gas backgrounds. The first three layers can be moved close to the beam, to maximize acceptance for heavy photons. A cooling system removes heat from the electronics, and cools the sensors to improve their radiation hardness. The sensors are readout with onboard sensors which pass signals to Front End boards for digitization and transmission out of the vacuum enclosure.

### 5.1.1 Hazards

Hazards to personnel include the high voltage which biases the sensors, and the low current which powers the readout electronics.

Hazards to the SVT itself include mechanical damage, radiation damage, and overheating. Hazards to the vacuum system could arise from excessive SVT outgassing or coolant leaks. SVT mechanical damage could occur if the top sensors are accidentally driven into the bottom sensors.

Radiation damage could occur in the SVT if the sensors are driven too close to the beam, the beam moves into the sensors, the beam interacts upstream to produce excessive radiation, or excessive beam currents create more radiation than can be tolerated.

Overheating can occur in the SVT if the cooling system is performing inadequately or if a cooling system leak develops.

### 5.1.2 Mitigations

Hazards to personnel are mitigated by turning off HV and LV power before disconnecting cables or working on the sensors and internal electronics. Hazards to the Hall B vacuum system have been mitigated by extensive testing of all components to ensure low-outgassing rates, construction of the SVT and electronics in a clean room, and tests



of the cooling system to high pressures to prove leak-tightness. The coolant used, water glycol, would not cause irreparable damage to the vacuum system if a leak occurred. If the system pressure increases, the coolant supply is halted.

Possible mechanical damage has been mitigated by designing the channels which hold the sensors to touch before any modules would touch. Software limits and limit switches on the motion controllers also prevent the sensors from moving into each other and too close to the beam.

Radiation damage from the beams is mitigated in several steps. First, beam size and halo must conform to beam requirements before beams are passed through the detector. Second, the beams are centered between the top and bottom sections of the SVT. Third, an upstream collimator is aligned with the "centered" beam position to intercept the beam if it moves off nominal position. Fourth, beam halo monitors sense a rise in backgrounds if the beam moves off nominal position, activating the FSD and removing the power permissive to the SVT. Fifth, precision movers position the SVT layers precise and safe distances from the beam. Finally, beam currents and target thicknesses are carefully chosen to avoid over-radiating the silicon sensors.

Overheating is mitigated by requiring good coolant flow, proper coolant temperature, good vacuum (assuring no coolant leakage), and sensor temperature in range in the interlock for SVT HV and LV power.

### 5.1.3 Responsible Personnel

Individuals responsible for the system are:

Name	Dept.	Phone	email	Comments
Tim Nelson	SLAC			First contact
Omar Moreno	UCSC			Contact
Sho Uemura	SLAC			Contact
Per Hansson	SLAC			Contact

Table 5.2: Personnel responsible for the silicon tracker.

## 5.2 Electromagnetic Calorimeter

The Electromagnetic Calorimeter (ECal) consists of 442 lead-tungstate ( $\text{PbWO}_4$ ) crystals with avalanche photodiode (APD) readout and amplifiers enclosed inside a temperature controlled enclosure. There are two identical ECal modules positioned above and below of the beam plane. In order to operate the calorimeter modules, high voltage and low voltage are supplied to each channel. The high voltage is  $< 450$  V and  $< 1$  mA. The required low voltage is  $\pm 5$  V for preamplifier boards. Constant temperature inside the enclosure is kept by running a coolant through the copper pipes that are integrated into the enclosure using laboratory chiller. Cooling system should provide temperature stability at the level of  $1^\circ\text{C}$ .

### 5.2.1 Hazards

Hazards associated with this device are electrical shock or damage to the APDs if the enclosure is opened with the HV on. There is also hazard associated with coolant leak that may damage preamplifier boards.

### 5.2.2 Mitigations

Whenever any work has to be done on the calorimeter, independent it will be opened or not, always turn HV and LV off. Turn chiller off if enclosure will be opened for maintenance. Any large (more than couple of degrees in C) must be investigated to make sure that there are no leaks.

### 5.2.3 Responsible Personnel

The authorized personnel is shown in table 5.3.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Stepan Stepanyan	Hall-B	7196		<a href="mailto:stepanya@jlab.org">stepanya@jlab.org</a>	<i>1st Contact</i>
Raphael Dupre	ORSAY	7252		<a href="mailto:dupre@ipno.in2p3.fr">dupre@ipno.in2p3.fr</a>	<i>2nd Contact</i>

Table 5.3: HPS calorimeter: authorized personnel.

## 5.3 CLAS12 Construction work

The main 12 GeV activity in Hall-B during the HPS engineering run will be assembly of the CLAS12 Torus magnet. Beam running would occur during evenings, nights, and weekends or during other periods when it would not conflict with the regularly scheduled assembly of the CLAS12 Torus coils.

### 5.3.1 Hazards

There are no personal hazards associated with the running over evenings, nights, and weekends, and continue with torus assembly during the normal work hours. The only hazard is possible delays of start of the torus work after beam the delivery due to possible activation of beamline parts close to the torus assembly fixtures.

### 5.3.2 Mitigations

Every time Hall will switch from beam running to torus assembly, a fully survey of the hall will be conducted and Hall will be brought to "Restricted Access". Normally this will happen very early in the morning of work day (6am). If elevated radiation near

torus assembly fixtures are found, work on torus must be delayed until conditions are acceptable.

However, we do not expect any excess radiation in the Hall or activation of any beam line components near the torus assembly area. The HPS target is located  $\sim 20$  meters upstream of the assembly area and only tuned electron beam, couple of hundred micron wide, will be transported in vacuum through hall to the target. If beam conditions are not acceptable, which may result excess radiation, beam tune will be performed. Every time beam tune is required the Hall-B tagger magnet will be energized and beam will be dumped on the Hall-B tagger dump, shielded hole in the floor  $\sim 15$  meters upstream of the torus assembly area.

### 5.3.3 Responsible Personnel

Individuals responsible for the coordination of the torus assembly and HPS engineering run:

Name	Dept.	Phone	email	Comments
PDL	Hall-B			1nd contact
S. Stepanyan	Hall-B		stepanya@jlab.org	contact
D. Kashy	Hall-B		Kashy@jlab.org	contact

Table 5.4: Personnel responsible for coordination of the HPS run and the torus assembly.

# Part III

## Slow Controls

# Chapter 6

## Overview <sup>1</sup>

The basic components of the system are:

- Input/Output Controllers (IOCs) - VME systems
- Operator Interfaces (OPI) - Computers capable of executing EPICS tools to interact with the IOCs. Some of the most used tools in Hall B are StripTools allowing to monitor the behavior of one or more signals as a function of time.

### 6.1 System's Components

The tasks assigned to the various IOCs are,  
ioch-1 Beam current/position monitors.  
ioch-2 Beam helicity asymmetry.....

### 6.2 Operating Procedures

Log into the Hall B control system through one of the computers....  
To start any of these applications, use .....

### 6.3 AlarmHandler

The “AlarmHandler” notifies the user when either a signal being monitored is outside some pre-defined limits or communication with the IOC in which the signal resides has been lost. “AlarmHandler” will only detect an abnormal signal condition if the signal is included in the application configuration file and, the corresponding IOC database record is set to produce an alarm condition.

---

<sup>1</sup>Authors: K. Livingston [kliv@jlab.org](mailto:kliv@jlab.org)

## 6.4 StripTool

Strip Tool plots a real-time strip chart of the values of one or more signals. It is useful to monitor data trends. A detail description of the options and operation of this application can be found in the Strip Tool Users Guide<sup>2</sup> with one difference; the version used by Hall A does not have a “print” function. To print a strip chart use the application “Snapshot” described below.

## 6.5 Snapshot

Snapshot refers to a KDE desktop application (ksnapshot) which allows to grab an image of either the whole screen or an individual window. The image can then be sent to a printer or stored on disk.

## 6.6 Troubleshooting Procedures

The status of most IOC's can be seen by opening the ....

Rebooting of the IOC's is accomplished in several ways depending on the specific IOC.

### 6.6.1 Authorized Personnel

The authorized personnel is shown in table 6.1.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Ken Livingston	Hall-B/Glasgow	XXX		<a href="mailto:kliv@jlab.org">kliv@jlab.org</a>	

Table 6.1: Slow controls: authorized personnel.

---

<sup>2</sup><http://www.aps.anl.gov/epics/extensions/StripTool/index.php>

# Part IV

## Data Acquisition and Trigger

# Chapter 7

## Spectrometer Data Acquisition <sup>1</sup>

The Hall B data acquisition uses CODA [6] (CEBAF Online Data Acquisition), a toolkit developed at Jefferson Lab by the Data Acquisition Group. Up to date information about the Hall B DAQ is kept at <sup>2</sup>.

We typically run with .....

The trigger supervisor is a custom-made module built by the data acquisition group. Its functions are to synchronize the readout crates, to administer the deadtime logic of the entire system, and to prescale various trigger inputs. We have a trigger supervisor...

The trigger management software is described in the Trigger chapter.

### 7.1 General Computer Information

In the counting room we have various computers for DAQ, analysis, and controls. The controls subnet is the responsibility of .....

### 7.2 DAQ checklist

Things to check before experiment starts

1. ....
2. raster cabling
3. check portservers
4. check reset
5. check trigger latch is connected
6. helicity if needed

---

<sup>1</sup>Authors: S.Boyarinov [boiarino@jlab.org](mailto:boiarino@jlab.org)

<sup>2</sup>[http://clasweb.jlab.org/equipment/daq/daq\\_trig.html](http://clasweb.jlab.org/equipment/daq/daq_trig.html)



- 7. scalers
- 8. synchronization time stamp

## 7.3 Beginning of Experiment Checkout

This section describes the checkout of DAQ and trigger needed before an experiment can start.

1. First ensure that all the fastbus, VME, CAMAC, and NIM crates are powered on.
2. Make sure the HV is on for all detectors and that the values are normal.
3. Start the scalers display following the instructions below and check that the rates from detectors are normal.
4. Startup CODA using the directions below and start a run. With the trigger downloaded and the HV on, you are taking cosmos data, typically at a rate of 3 Hz per spectrometer.

## 7.4 Running CODA

This section describes how to run CODA for the CLAS12 DAQ.

Here is how to start and stop a run. Normally, when you come on shift, runcontrol will be running. If not, see the section on “Cold Start” below. To start and stop runs, push the buttons “Start Run” and “End Run” in the runcontrol GUI. To change configurations use the “Run Type” button. If you have been running you will first have to push the “Abort” button before you can change the run type. Typically the configurations you want are the following.

A note about pedestal runs. They have the exclusive purpose of obtaining pedestals used for pedestal suppression. For details about what is done and hints for getting pedestals for analysis (which does not want the PEDRUN result), see [/ped/README](#).

### 7.4.1 Some Frequently Asked Questions about DAQ

- *Q: Where is the data ?* Use a command ....

Files are archived automatically to tape in the MSS tape silo. Two tape copies are made. Data are purged from disk automatically. Users should *never* attempt to copy, move, or erase data.

- *Q: How to adjust prescale factors ?* Edit the file .....
- *Q: What is the deadtime ?* The deadtime is displayed in ....

- *Q: Where are the crates ?* Fastbus crates ROCXXX are ...
- *Q: Why is the deadtime so high ? (and related)* Search for answers among the following....

### 7.4.2 Quick Resets

Problems with CODA can usually be solved with a simple reset. If not, try a Cold Start (see next section). Do not waste an hour of beam time on resets; if they fail, call an expert. The expert claims he can restart CODA 90% of the time within 10 minutes.

If a ROC (ReadOut Controller, or crate) is hung up, reboot by going the workspace “Components” and typing ..... If this doesn’t work, try pressing the reset button which is on the “Crate Resets” section of the

### 7.4.3 Cold Start of CODA

If CODA is not running, or if it gets hung up, you can do a cold start. Frequently a subset of these steps is sufficient to recover from a hangup, but it takes some experience to realize the minimum of steps that are necessary, so the simplest thing is to do them all, which takes a few minutes.

- Make sure the fastbus and VME crates are running. The crates are usually known by .....
- Start runcontrol and the other necessary processes by typing .....
- In runcontrol, press the “Connect” button. Wait 5 seconds and press “Run Types”. After configure and before download, press the “Reset” button in the upper left corner. Choose the run type from the dialog box (see section on Running CODA for descriptions of run types).
- After you configure and download the Run Type, you can “Start Run” to start a new run.

### 7.4.4 Recovering from a Reboot of Workstation

If the workstation from which you are running CODA was rebooted, here is how to recover DAQ. Login as the relevant account, which is usually a-onl for 1-DAQ operation. Passwords for the online accounts should be available on a paper on the wall in the counting room, or ask the run coordinator. In the workspace for “Components” telnet into all the ROCs. If the x-terms windows are not available, type .....

## 7.5 Electronic Logbook and Beam Accounting

Two tools are available for logging information by the shift workers: (1) The Electronic Logbook, and (2) The Hall Beam–Time Accounting Table.

The electronic logbook is a web-based repository of logbook data.

The Hall Beam–Time Accounting Table is the mechanism to summarize and record how the beam time in a shift was spent. The shift leader is responsible for submitting this table at the end of the shift. When submitted, the data are logged in a database and a summary is e-mailed to various people like the run coordinators and the hall leader.

## 7.6 Authorized Personnel

The authorized personnel is shown in table 7.1.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Sergei Boiarinov	Hall-B	5795		<a href="mailto:boiarino@jlab.org">boiarino@jlab.org</a>	<i>Contact, first on call</i>

Table 7.1: DAQ: authorized personnel.

server IP	Port	Device
-----------	------	--------

Table 7.2: Port Servers for DAQ

# Chapter 8

## Trigger Hardware and Software

1

### 8.1 Overview

Here we give a brief overview of the hall A trigger, including its hardware arrangement, the logic of the trigger, and the usage of the software control. Diagrams of the hardware layout are shown in accompanying figures.

The trigger design is quite flexible and it is relatively easy to add detectors to define new trigger types or to modify existing ones, so long as the detector is fast enough. The trigger supervisor also allows for the possibility of 2nd level triggers which could be used for a later decision.

### 8.2 Components

The trigger schematics is shown in Fig. 8.1.

Here we describe the software control of the modules involved in the trigger.

Here are the instructions to download the trigger.

The user should look for suspicious error messages in the window from which trigsetup was launched, e.g. to check if connection to the crate is ok.

If individual modules need to be modified for test purposes etc. (e.g. to change thresholds), one may use the expert mode.

#### 8.2.1 Authorized Personnel

The authorized personnel is shown in table 8.1.

---

<sup>1</sup>Authors: S. Boiarinov [boiarino@jlab.org](mailto:boiarino@jlab.org)

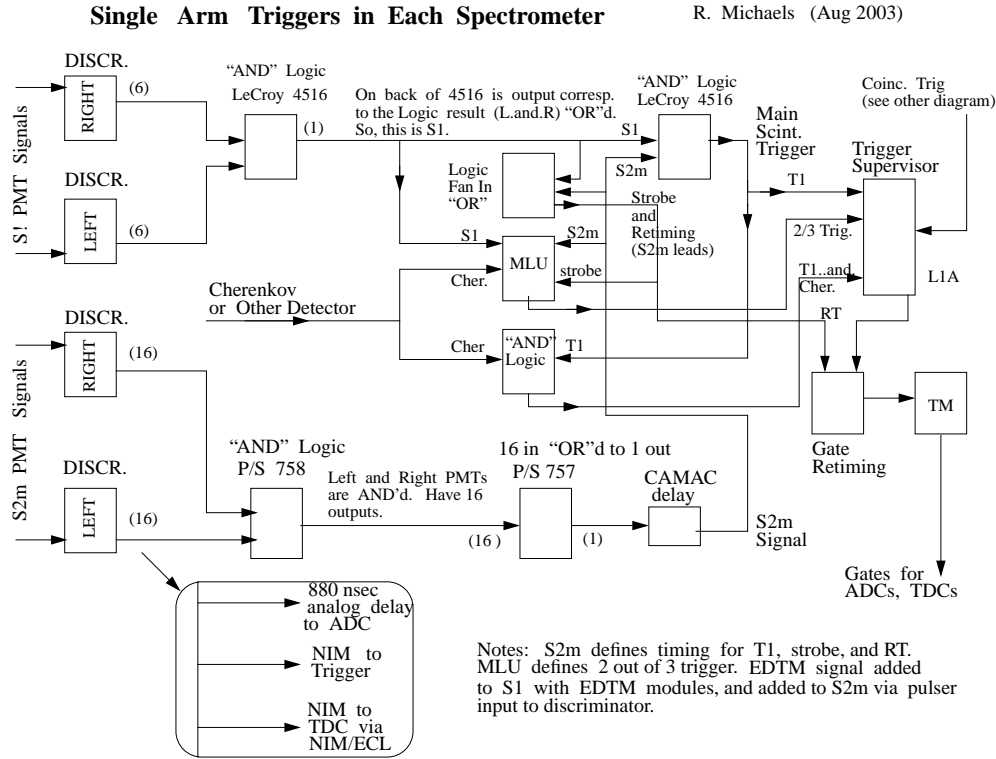


Figure 8.1: Trigger Circuit.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Sergei Boiarinov	Hall-B	5795		<a href="mailto:boiarino@jlab.org">boiarino@jlab.org</a>	<i>Contact</i>
ValeryKubarovsky	Hall-B	5647		<a href="mailto:vpk@jlab.org">vpk@jlab.org</a>	

Table 8.1: Trigger: authorized personnel.

# Chapter 9

## Online Analysis, Data Checks

<sup>1</sup>

The following tools are available for checking data online.

### 9.0.2 Scaler Display and Scaler Events

Scaler rates and values are displayed using a ROOT based display.

Normally this is already running, but if it is not running, login as .... and go to the appropriate directory by typing .....

## 9.1 Analysis using CLAS12 Monitoring

CLAS122 offline software package for analyzing Hall B data. This code is documented in a separate chapter but it is worth mentioning here in a list of essential tools for checking data online.

### 9.1.1 Responsible Personnel

The responsible personnel is shown in table 9.1.

Name (first,last)	Dept.	Call [5]		e-mail	Comment
		Tel	Pager		
Sergei Boiarinov	Hall-B	5795		<a href="mailto:boiarino@jlab.org">boiarino@jlab.org</a>	<i>Contact</i>
Veronique Ziegler	Hall-B	6003		<a href="mailto:ziegler@jlab.org">ziegler@jlab.org</a>	

Table 9.1: Online analysis: authorized personnel.

---

<sup>1</sup>Authors: S.Boiarinov [boiarino@jlab.org](mailto:boiarino@jlab.org)

# Part V

## Offline Analysis Software

# Chapter 10

## CLAS12 Offline Analysis

<sup>1</sup>

The standard offline analysis software for Hall B data is ....

The XXX code performs tracking in the focal plane and reconstruction to the target. The tracking algorithm has been shown to be accurate for events with one cluster per plane. Noisy events with higher cluster multiplicity and events with more than one good track in the focal plane may not be reconstructed correctly by the present version of the code, but work is in progress to make this type of analysis also reliable.

The scintillator, Cherenkov, and shower counters classes perform basic decoding, calibration (offset/pedestal subtraction, gain multiplication), and summing (for Cherenkovs) or cluster-finding (for showers) of hits. The cluster-finding algorithm of the shower class is basic and currently only capable of finding a single cluster per event. These classes are largely generic and should be able to accommodate most new detectors of the respective type, even with a different geometry and number of channels.

Data of interest can be histogrammed and/or written to a ROOT Tree in the output file. The contents of the output is defined dynamically at the beginning of the analysis. Both 1- and 2-dimensional histograms are supported. Histograms can be filled selectively using logical expressions (cuts).

Table XXX lists the analysis modules available

---

<sup>1</sup>Author: V. Ziegler [ziegler@jlab.org](mailto:ziegler@jlab.org)



# Bibliography

- [1] EPICS Documentation. WWW page. URL <http://www.epics.org/>.  
see also <http://www.aps.anl.gov/asd/controls/epics/EpicsDocumentation/WWPages/EpicsDoc.html>.
- [2] JLab. *EH&S Manual*. URL <http://www.jlab.org/ehs/ehsmanual/>.
- [3] JLab. *Personnel Safety System (PSS) manual*. URL [http://www.jlab.org/accel/ssg/user\\_info.html](http://www.jlab.org/accel/ssg/user_info.html).
- [4] *Accelerator Operations Directive*. URL [http://opsntsrv.acc.jlab.org/ops\\_docs/online\\_document\\_files/ACC\\_online\\_files/accel\\_ops\\_directives.pdf](http://opsntsrv.acc.jlab.org/ops_docs/online_document_files/ACC_online_files/accel_ops_directives.pdf). *URL is available inside JLab site.*
- [5] Jefferson Lab, (12000 Jefferson avenue, Newport News, VA 23606). URL <http://www.jlab.org>. Telephone numbers: (757)-269-XXXX, Pager numbers: (757)-584-XXXX.
- [6] CODA documentation. WWW page. URL <http://coda.jlab.org/>.